USER MANUAL / MANUAL DE USUARIO / MANUAL DE USUÁRIO / MANUEL DE  
L'UTILISATEUR / MANUALE DELL'UTENTE / BENUTZERHANDBUCH / GEBUIKREK-  
SHANDLEIDING - BAS/AU1536

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Medición de la  
masa ósea

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AURA

# Báscula inteligente digital **Wireless BT**

## Check your phone model

Please confirm your phone or PC tablet can be adapted before using APP

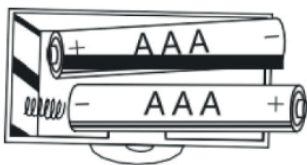
Bluetooth 4.0 + and iOS 7.0 +

Bluetooth 4.0 + and Android 4.3 +

## Installing the battery

1. Remove the battery cover at the back of scale
2. Insert 2 X 1.5V AAA batteries (Included) into the scale correctly
3. Put back the battery cover.

## Set the unit

To set “unit”, press the button located on the back of scale, right behind the display. Choose the right unit you need: Kg, Lb or St.

## Installing APP

- Switch on Bluetooth
- Phone or tablet connects to Net (WiFi, 4G or 3G)
- Use your device to scan QR code or search “OKOK international” from App Store or Google Play to install App

iOS

Android

If Android can't download the "OKOK international" directly through scan barcode, you can get into "Google play" to download it.

## **Use the devices and register, complete user informations**

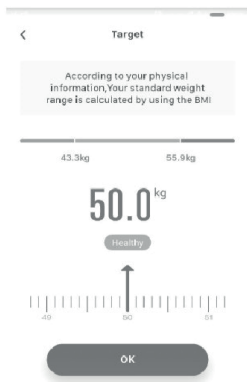
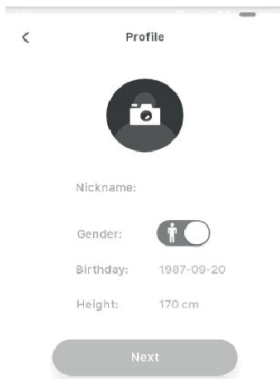
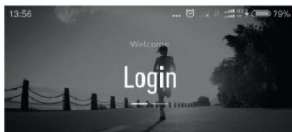
Open "OKOK International"

1. If you have the account, fill the email and password to enter;
2. Haven't the account, click "Join us", register your personal information correctly.

Nickname > Gender > Birthday > Height > Target

(Target: According to your physical information, your standard weight range is calculated by using the BMI)

3. If you are registered and forgot the Password, select "Forget password".

Set the phone open "Bluetooth"  
Enter into the main screen on  
APP

Click the “Device”; Turn on the scale; Bind the Bluetooth scale on APP

Device

Bind Bluetooth Scale

Bluetooth not turned on, please turn it on

Please step scale lightly and  
keep its screen light on.

Not Bind Now

Click the “Go” add the device.

Device

Add device

# Measuring

Keep the Bluetooth on and stand on the scale, the value will display on your screen

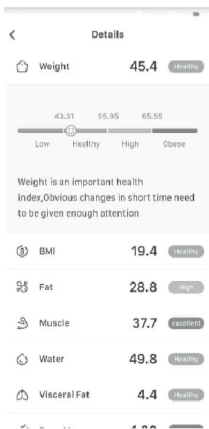
07/25/2018 12:38

History

Weight	50.0	Healthy
BMI	21.4	Healthy

Click the recorded value, view measurement results

Click the "Trend"  
Easy to check each item data.

## APP Functions

Click the image icon to go to the personal profile, you can change / monitor the personal information accordingly.

Displays weight, body fat or BMI, you can scroll down to see histories.

Click “weight” area to see your Weight, Fat (%), Muscle (%), Hydration (%), BMI, Visceral Fat (%), Bone (LB), Calories (Kcal/day).

Your personal information is related with Algorithm accuracy, thus please make sure that your personal details like gender, height, age are entered correctly.

Note:

1. If only weight data is available, click “weight” will simply pop up weight and BMI data.
2. The report is based on full evaluation of all measurements thus it will function only when all measurements are valid, and all data are automatically saved while weighing.
3. If your phone is not available for pairing when you are weighing, the data will be updated when they connect.

It will display you maximum, minimum, and average weight, body

fat percentage, hydration, muscle mass percentage, bone mass percentage in week, month, quarter periods.

The generated graphs help you easily see the trend of this data.

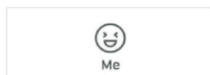
Open the “OKOK” APP, click “Go” to add the device before to weight, otherwise can't display the data same as scale.

1. Open the bluetooth function on phone, turn on the scale, keep its screen light on, Bind the Bluetooth scale

2. If one scale have connected the “OKOK” APP by phone, now want to change another phone use APP to connect the scale, firstly, Unbind the bluetooth scale.

1.  **My target weight** Editor the target weight here.
2.  **Family Member**

Invite up to other 7 household members add his / her name, gender, date of birth & height.

3.  Setting interface used to set account and

unpair the scale and so on

## Main Function

- Display: All data displayed when power on
- Capacity: 6 Kg - 180 Kgs
- Division: 1D = 100 g / 0.2 lb
- Unit: KG / LB / ST
- Voltage: 3.0 V
- Battery: 2 X AAA battery
- Best working temperature: 5 ~ 35° C
- Storage temperature -10 ~ 60° C
- Moisture environment: < = 90% RH
- Bluetooth frequency band: 2.4 Ghz
- Maximum emitted RF power < 100 mW

## Special Display

- Low battery less than 2.4V, display “Lo” and power off
- Overload indication: overload more than 2.5 Kgs, total 182.5 Kg display “Err”, fat and other data error display “Err”
- The auto power off: no operation about 20 seconds will automatically be turned off

## Special Display

1. Keep scale surface dry to avoid accident by slipping.
2. Please put the scale in a dry place to protect the internal electronic components.

3. Please use micro cloth to clean the scale, not immersed in water.
4. Please don't remove the scale except replace the battery.
5. Please don't use heavy percussion or falls from a height.
6. Please don't persistent put heavy on the scale to avoid damage to the high precision sensor.

## Main Function

1. Please keep the scale away from wet floor, don't use the scale
2. Please stand on the middle of scale in case of turnover.
3. The platform of the scale is made of glass, please avoid sharp things scratching the surface, excessive forces or sudden shocks, unexpected fallen-off will break the glass and injure the user.
4. The Scale is not waterproof. Keep the surrounding environment dry, do not submerge it in water.
5. use cloth with neutral cleaning liquids to clena the surface, keep away contacting with organic solvent nor chemicals.
6. When the batteries are in low power status, please remove the batteries in case for leaking liquids corrode the scale.

## Notices

1. Please stand on the scale in the right direction.
2. Please stand on the scale no less than 10 seconds if you measure body fat.
3. Please confirm gender, height, weight data is correct.

## GUARANTEE

### DEVICE

ITEM: \_\_\_\_\_

Serial number: \_\_\_\_\_

### ESTABLISHMENT

Purchase date: \_\_\_\_\_

Establishment signature

+34 91 713 94 98

**INNOVA CELULAR SL**  
C/ Severo Ochoa N° 9 nave 8B  
28521 Rivas-Vaciamadrid  
Madrid

**Note:** Do not expose the scale at extreme temperatures, the batteries could exploit.

**NOTE:** The BASCULA INTELIGENTE DIGITAL WIRELESS BT (BAS/AU1536 (BLUETOOTH SCALE EBS-016, EBS-016K, EBS-016W, EBS016R, EBS016G, EBS016A, EBS016B, EBS016C) has a three-year warranty from the date of purchase, except for the batteries, which must be maintained by the customer.

## Declaration of conformity

Through this INNOVA CELLULAR, it declares that the BASCULA INTELIGENTE DIGITAL WIRELESS BT (BAS/AU1536 (BLUETOOTH SCALE EBS-016, EBS-016K, EBS-016W, EBS016R, EBS016G, EBS016A, EBS016B, EBS016C) complies with the essential requirements and any other applicable or enforceable provisions of Directives 2014/53/EU (RED), RoHS 2011/65/EU, annex II (EU) 2015/863 and mentioned en (EU) 2017/2102 , of the European Parliament and of the Council, transposed to Spanish legislation through Royal Decree 188/2016, of May 6, 2016. To see the complete declaration of conformity, you can access it through the following link:

[http://www.innovacelular.com/images/web/descargas/descargas/declaracion/declaracion\\_bas-au1536.pdf](http://www.innovacelular.com/images/web/descargas/descargas/declaracion/declaracion_bas-au1536.pdf)

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**Working hours:** Mon to Thu from 08:00 to 14:00 and from 15:30 to 17:30 and Fri. from 09:00 to 14:00

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